

AIR~SUPPORT RIBCUSHION

Alternating air mattress system



The Hawksley Ribcushion mattress system is an effective and flexible means of reducing pressure ulcers and sores. The Air~Support mattress is priced at a level that makes it economic to be put on every long stay bed and will deliver huge savings by reducing pressure injury risks for your patients.

- Ribbed alternating two-chamber pressure mattress
- Laser-drilled micro pores for controlled air loss and patient comfort
- Tough, comfortable, easy-to-clean material
- Reliable, quiet and compact air pump system
- Suitable to be hung on a standard NHS bed
- Optional CPR valve available
- Capable of supporting patients up to 150kg
- Effective for reduction of pressure ulcers, venous insufficiency and swollen limb conditions
- Rental and lease options available

AIR~SUPPORT RIBCUSHION

Alternating air mattress system

Technical specifications

Ribbed mattress

Colour	Green with green cover
Height inflated	120mm (4.5")
Number of ribs	18 ribbed segments
Dimensions	1981mm L x 914mm W (78" X 36") c/w tuck-in flaps
Material	Non-toxic stain-resistant nylon with PVC, 0.41mm thick
Warranty	12 months

Air pump

Cycle time	6 minutes
Pressure range	40mm Hg to 130mm Hg
Power cord	1.8 metres
Colour	Light grey/green base
Dimensions	11cm height x 22.5cm deep x 10cm wide
Weight	1.3kg
Construction	Impact-resistant plastic shell and base
Warranty	12 months


Ordering information

39000-00	Air~Support Ribcushion air bed 220V c/w mattress
39001-00	Ribcushion mattress only




Alternating pressure cycles

Mattress cycles by inflating the A and B cells shown below in a total cycle time of 6 minutes broken into four phases.

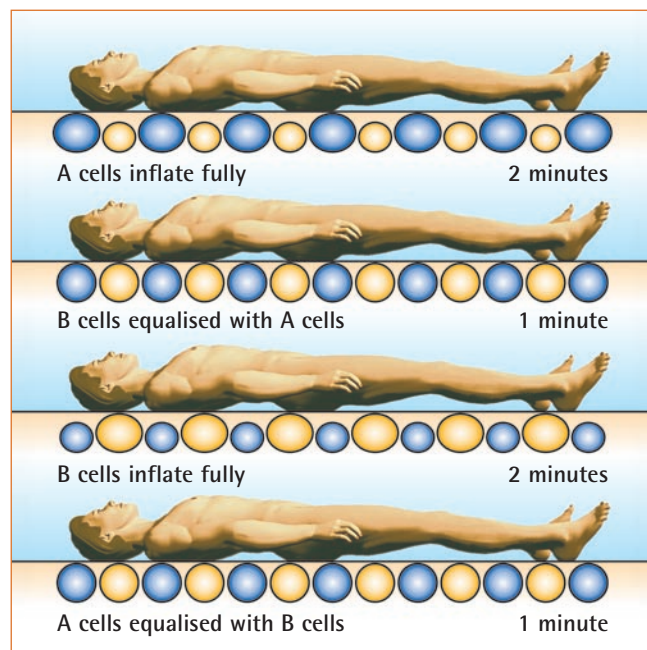

Typical weight distribution of adult male

Alternating pressure relief redistributes the pressure on the patient to allow better blood circulation.

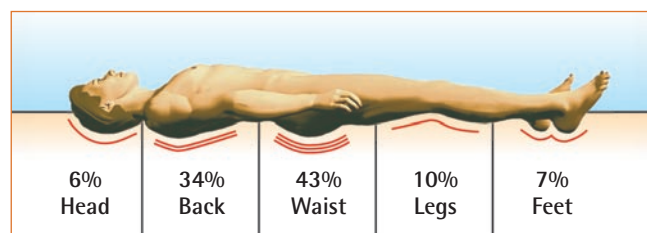

Hawksley • Medical and Laboratory Equipment

Marlborough Road, Lancing Business Park, Lancing, W. Sussex BN15 8TN, UK

T +44 (0)1903 752815 • F +44 (0)1903 766050 • www.hawksley.co.uk

